Installationsanleitung
Installation manual

KB/P511

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Reagent cooler for TresCon Pges

For your safety


Caution
For safety reasons devices standing on the reagent tray of TresCon must not be operated with line power (fluids!).

- **Operate the cooler using the enclosed 12 V power supply only.**
- **Install the power supply protected from fluids and apart from the reagent tray.**

Scope of delivery


Component	Quantity
Cooler, 12 VDC, with connection cable	1
Power supply unit (input 220-240 V / 50 Hz, output12 VDC) with line cord	1
Tube bushings, preassembled, with assembly material	4

Assembly

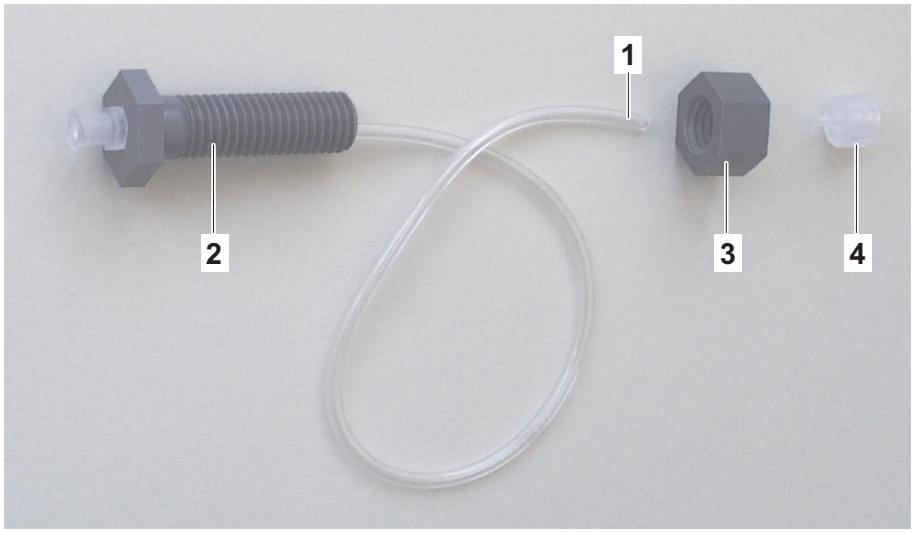1	Mount the four tube bushings in the borings of the lid: – Take off the lid. The lid can be taken out of the hinge for easier handling. – Insert the tube ends (1) into the borings from the outside and insert the bushings (2) into the borings. – Fix the bushings on the inside using the nuts (3). – Plug a Luer connector (4) on each free tube end.
2	Place the reagent containers A to D in the cooler under the bushings that are marked correspondingly and connect the tube ends with the screw caps.

Lid with tube bushings screwed in

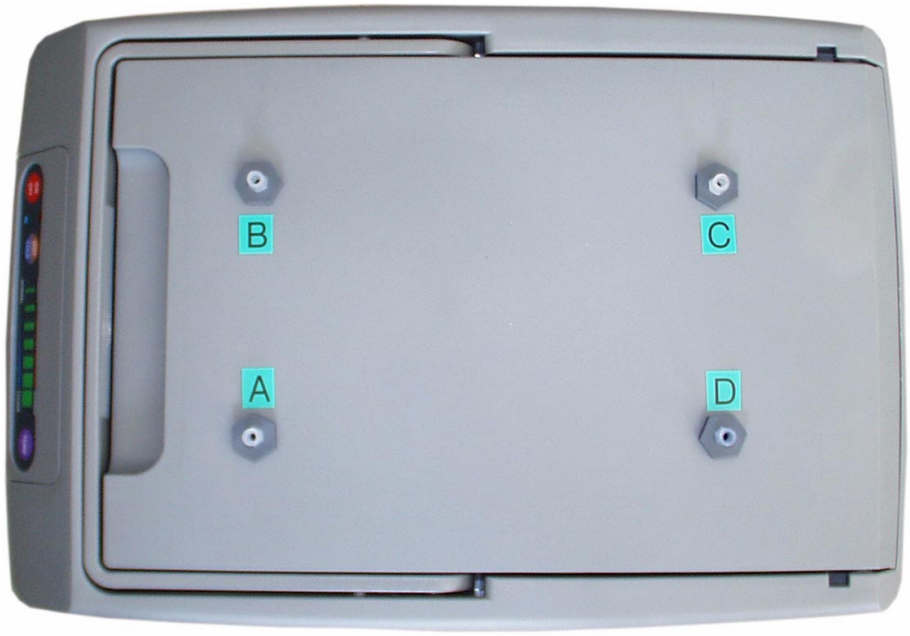cooler with assembled tubes


3	Close the lid and connect the TresCon reagent tubes A to D with the correspondingly labeled bushings.

Electrical connection

4	Plug the power plug of the power supply unit into a power socket. Make sure the values for line voltage and line frequency stated on the nameplate agree with the values of the power socket.
5	Install the power supply protected from fluids and apart from the reagent tray.
6	Using the connection cable, connect the cooler to the cigartette lighter socket of the power supply unit. The cooler is now ready to operate.

Commissioning

7	To switch on the cooler, press the ON/OFF key.
8	Adjust the operation of the cooler so that the inside temperature is approx. 10 °C. To do so, select the operating mode, "cooling", by actuating the "HOT/COLD" key (blue LED is illuminated). Subsequently press the "TEMP" key repeatedly until altogether 5 green LEDs are illuminated on the "NORMAL < COLDER" status display.

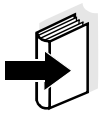
Note
For more information on the operation of the cooler please refer to its operating manual.

What to do if ...

Cooler does not work (no display, ventilator does not rotate)

Cause	Remedy
Power supply interrupted	– Check connections – Check the fuse of the power supply unit and replace it if necessary – Check the fuse in the connection cable of the cooler and replace it if necessary. The fuse is in the connector plug (cigarette lighter plug).

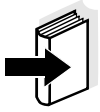Note
For more information on error elimination please refer to the operating manual of the cooler.