QUICK START GUIDE

ba77019y02 03/2014

pH/ION 3310

WTW

a xylem brand

Copyright © 2016 Xylem Analytics Germany GmbH
Printed in Germany.

pH/ION 3310 - Contents

EN English 4
DE Deutsch 6
CS Česky 8
DA Dansk 10
EL Ελληνικά 12
ES Español 14
FI Suomi 16
FR Français 18
HU Magyar 20
IT Italiano 22
JA 日本語 24
KO 한국어 26
NL Nederlands 28
NO Norsk 30
PL Polski 32
PT Português 34
RO Română 36
RU Русский 38
SK Slovenský 40
SV Svenska 42
TR Türkçe 44
ZH 中文 46

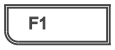
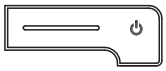
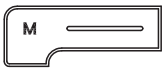
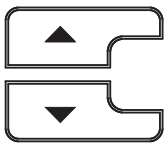
pH/ION 3310

Overview


The pH/ION 3310 meter enables you to perform measurements (pH) quickly and reliably.
The pH/ION 3310 meter provides a maximum degree of operating comfort, reliability and measuring certainty.

Keypad

Keys are indicated by brackets <..>.
The key symbol (e.g. **<ENTER>**) generally indicates a short keystroke (under 2 sec). A long keystroke (approx. 2 sec) is indicated by the underscore behind the key symbol (e.g. **<ENTER__>**).


  **<F1>:** **<F1__>:** **<F2>:** **<F2__>:**	Situation dependent function (is shown on the display) e.g.: **<F1>**/[Menu]: Opens the menu for measurement settings **<F1__>**/[Menu]: Opens the menu for system settings
 **<On/Off>:**	Switches the meter on or off
 **<M>:**	Selects the measured parameter
 **<CAL>:** **<CAL__>:**	Calls up the calibration procedure Displays the calibration data
 **<STO>:** **<STO__>** :	Saves a measured value manually Opens the menu for the automatic save function
 **<RCL>:** **<RCL__>:**	Displays the manually stored measured values Displays the automatically stored measured values
 **<▲><▼>:**	Menu control, navigation
 **<ENTER>:** **<ENTER__>:**	Opens the menu for measurement settings / confirms entries Opens the menu for system settings
 **<AR>:**	Freezes the measured value (HOLD function) Switches the AutoRead measurement on or off

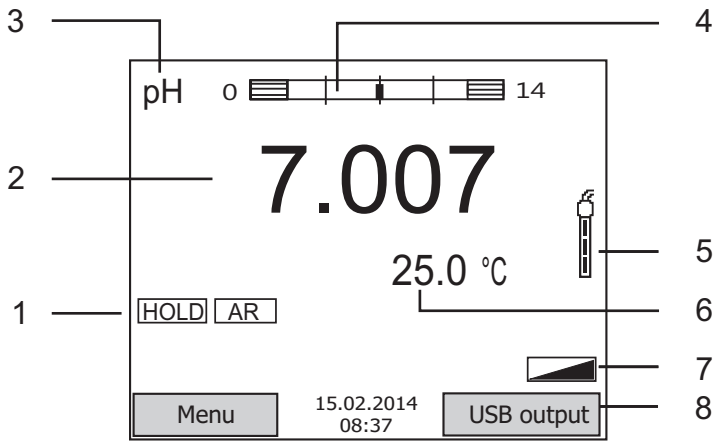Display

1	Status information
2	Measured value (with unit)
3	Measured parameter
4	Continuous measurement control (CMC function)
5	Sensor symbol (calibration evaluation)
6	Measured temperature (with unit)
7	Further status information
8	Softkeys and date/time

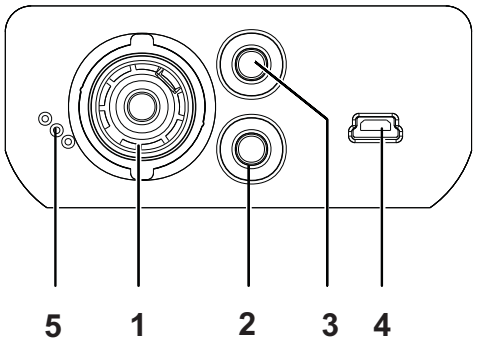Jack field /connections

1	pH/ION electrode
2	Reference electrode
3	Temperature sensor
4	USB B interface (device)
5	Service interface

Initial commissioning

1	Insert the supplied batteries. Make sure that the poles are in the correct position.
2	Switch on the meter with the **<On/Off>** key.

Caution

Follow the safety instructions of the used sensors.
The operating manuals of the sensors are available on the supplied CD and on the Internet under www.WTW.com.

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Xylem |ˈzīləm|

1) The tissue in plants that brings water upward from the roots;
2) a leading global water technology company.

We're a global team unified in a common purpose: creating advanced technology solutions to the world's water challenges. Developing new technologies that will improve the way water is used, conserved, and re-used in the future is central to our work. Our products and services move, treat, analyze, monitor and return water to the environment, in public utility, industrial, residential and commercial building services settings. Xylem also provides a leading portfolio of smart metering, network technologies and advanced analytics solutions for water, electric and gas utilities. In more than 150 countries, we have strong, long-standing relationships with customers who know us for our powerful combination of leading product brands and applications expertise with a strong focus on developing comprehensive, sustainable solutions.

For more information on how Xylem can help you, go to www.xylem.com.

Service and Returns:
Xylem Analytics Germany
Sales GmbH & Co. KG
WTW
Am Achalaich 11
82362 Weilheim
Germany

Tel.: +49 881 183-325
Fax: +49 881 183-414
E-Mail wtw.rma@xylem.com
Internet: www.xylemanalytics.com

Xylem Analytics Germany GmbH
Am Achalaich 11
82362 Weilheim
Germany