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Instruction manual for the
oxygene electrodes

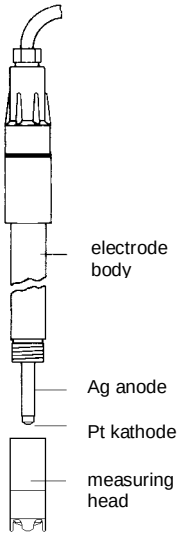
fig.1

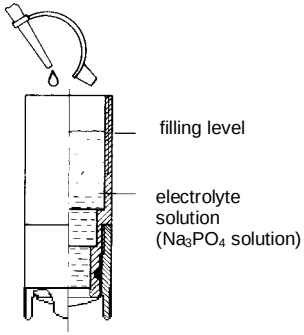
fig.2

Xylem Analytics Germany GmbH

Am Achalaich 11
82362 Weilheim
Germany

SI Analytics
Tel.: +49.(0)6131.66.5119
Fax: +49.(0)6131.66.5001
E-mail: Service-Sensors.si-analytics@xylem.com
www.XylemAnalytics.com

SI Analytics
a xylem brand

CE UK CA

Preparation

Attention!
With unsrewed measuring head the sensitive glass part at the end of the electrode and the PTFE membrane (0.025 mm) should not be destroyed. Both parts are good protected when assembled. An electrode with defective glass part must be replaced.

Filling of the measureing head with electrolyte solution

- Unscrew the measuring head
- Fill in electrolyte solution, about 30 drops (s. fig. 2)
- Remount the measuring head
- The electrode is ready to perform a measurement

Calibration

Please also note the instruction manual of the measuring device during measurement.

a) Zero point

- Fill the zero point solution from the glass ampoule into a beaker, connect the electrode to the oxygen meter and switch to setting „mg/l“ (polarisation).
- The zero point solution is free of oxygen, so after a few seconds the display reading is close to 0.0 mg/l.
- Turn the control NULLPUNKT/ZERO to read 0.0 mg/l

After replacement of the measuring head the first polarisation could last up to 15 min.

b) Slope

- Rinse the electrode with distilled water, expose it to air and wait for steady signal value (10 – 30 sec.).
- Read the temperature
- Adjust the control knob STEILHEIT/SLOPE to read the saturation value corresponding to the measuring temperature (see table below).

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Maintenance

If there is no more electrolyte in the measuring head, at first rinse the electrode with distilled water, then rinse the electrode with new electrolyte solution (s. fig. 2). In case of membrane defects the measuring head should be replaced by a new one.

Storage

If the electrode is not in use, the watering cap, filled with distilled water, should be mounted and the electrode stored in its original box.

Quality

Oxygen electrodes are manufactured with high precision by strict quality specifications and tested individually. Is there any fault observable on the originally packed electrode then the electrode must be returned immediately.

Supply

We supply each oxygen OX 1100 electrode with:

50 ml oxygen solution L 6708
Ampoule with zero point solution
Spare mounted measuring head

Spare parts

Besides the delivered service set you can order the following spare parts:

OX 920 50 ml oxygen electrolyte
OX 921 30 ml cleaning solution
OX 923 3 measuring heads
OX 925 maintenance kit (OX 920, OX 921, OX 923, SF 300)
OX 928 5 measuring heads, 50 ml oxygen electrolyte, 6 ampoules with zero point solution
OX 929 5 measuring heads
OX 060 60 ampoules zero point solution
L 6708 50 ml oxygen electrolyte
SF 300 abrasive sheet
Z 472 Watering cap

Further information can be found in our product catalogs.

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